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
ESTEREOTÁXICO AIMSISTEM TM03B

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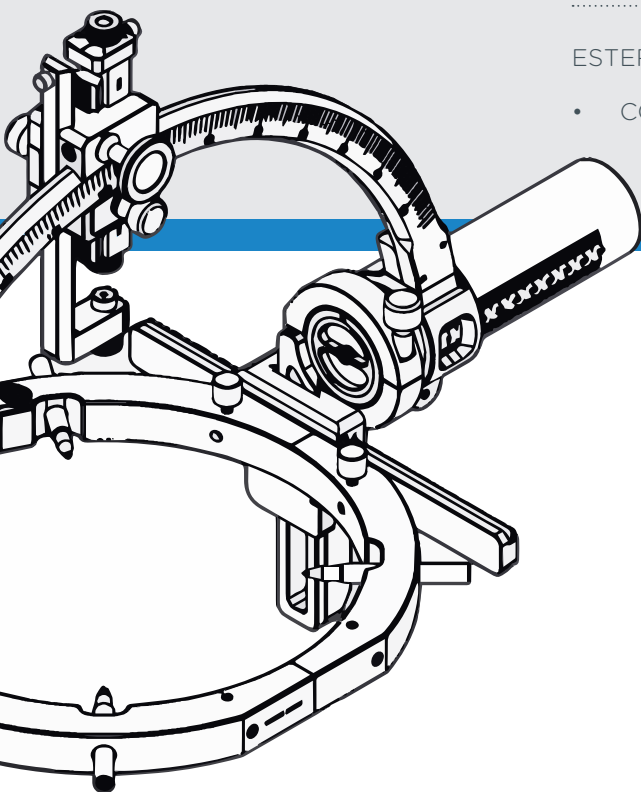M I C R O M A R

MANUAL DO USUÁRIO	3
USER MANUAL	45
MANUAL DEL USUARIO	87

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ESTEREOTÁXICO AIMS SYSTEM TM03B

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ÍNDICE

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STEREOTACTIC AIMS SYSTEM TM03B

- BASIC COMPONENTS

1. INTRODUCTION.....	48
2. ILLUSTRATIONS.....	48
3. COMPATIBILITY.....	48
4. USE OF EQUIPMENT	48
5. PRESENTATION	48
6. BASIC CONCEPTS AND OPERATION PRINCIPLE	49
7. AIMS SYSTEM TM03B.....	50
8. INTENDED USE.....	50
9. PROCEDURES	50
9.1. COMPONENT LIST	51
9.2. PREPARING THE PATIENT FOR THE SURGICAL PROCEDURE.....	53
9.2.1. DIRECT FIXATION.....	54
9.2.2. INDIRECT FIXATION.....	56
9.2.3. PREPARATION FOR IMAGING EXAM	60
9.2.4. MANUAL CALCULATION OF STEREOTACTIC COORDINATES.....	61
9.2.5. COMPONENT ASSEMBLY	69
10. SCALES.....	78
11. TRAYS FOR TRANSPORT / STORAGE / STERILIZATION.....	78
12. CLEANING AND STERILIZATION.....	79
13. LUBRICATION.....	80
14. PREVENTIVE AND CORRECTIVE MAINTENANCE	81
15. SYMBOLS	82
16. CARE AND PRECAUTIONS	83

PRODUCT DESCRIPTION:

STEREOTACTIC AIMSISTEM TM03B

MODEL: TM03B

ONLY USE THE EQUIPMENT AFTER CAREFULLY READING THE USER'S MANUAL AND RESPECTIVE WARNINGS.
CHECK IF THIS VERSION OF INSTRUCTION FOR USE IS ACCORDING TO THE PRODUCT PURCHASED.

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INFO@MICROMAR.COM

WWW.MICROMAR.COM

MICROMAR IND. COM. THE RIGHT TO MAKE CHANGES TO THIS MANUAL WITHOUT PRIOR NOTICE IS HEREBY RESERVED. CONTACT OUR SALES DEPARTMENT TO RECEIVE THE LATEST VERSION OF THE USER MANUAL.

THIS MANUAL CONTAINS INFORMATION PROTECTED BY COPYRIGHT

THIS MANUAL MAY NOT BE REPRODUCED, TRANSLATED OR TRANSMITTED BY ANY ELECTRONIC OR OTHER MEANS, IN WHOLE OR IN PART, WITHOUT EXPRESS AUTHORIZATION FROM MICROMAR

1. INTRODUCTION

This manual aims to provide user instructions on the safe and effective use the Stereotactic AimSystem TM03B. Before using the product (TM03B and its Accessories), the user should carefully read the manual and gain a total understanding of its contents, in addition to undergoing the training provided by Micromar. This manual should be kept with the equipment for reference in the event that questions arise.

2. ILLUSTRATIONS

The images in this manual illustrate the product and its accessories. There may be differences compared to prior equipment models and versions. These differences in illustrations do not interfere in its use.

3. COMPATIBILITY

The equipment described in this manual may only be used in combination with those components expressly recognized by Micromar as compatible. Contact Micromar for more information.

4. USE OF EQUIPMENT

The Stereotactic AimSystem TM03B is made to be assembled and used according solely to the procedures indicated in this manual and for the purpose for which it was developed.

5. PRESENTATION

The basic setup of the Micromar's Stereotactic AimSystem TM03B is shown, along with the components needed to perform the surgical procedures for which the product is developed.

The Stereotactic AimSystem TM03B has a line of accessories developed to aid the user during procedures, aimed at assuring surgical efficacy and patient safety.

6. BASIC CONCEPTS AND OPERATION PRINCIPLE

Etymologically, stereotaxy (from Greek, stereos - solid, three-dimensional; taxis - placement, position) means guided three-dimensional localization. Stereotactic neurosurgery means reaching a target on a region of the brain, using Cartesian coordinates that are defined in advance using a geometric orientation device, which is attached to the head and allows for an instrument to be guided to a target.

The Stereotactic AimSystem TM03B allows for targets to be located using Cartesian coordinates (X, Y and Z), as shown in Figure 1.

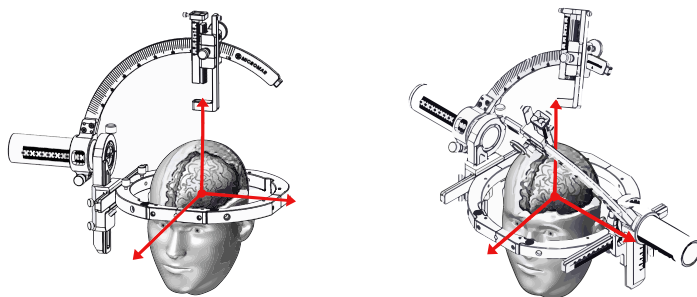
Figure 1

"X" = Right and left sides

"Y" = Back and front

"Z" = Top and bottom

The equipment uses an isocentric concept that allows a target to be reached from any point of entry; in other words, the target is reached with precision regardless of Alpha and Beta angles. Figure 2 shows how the equipment can be angled.

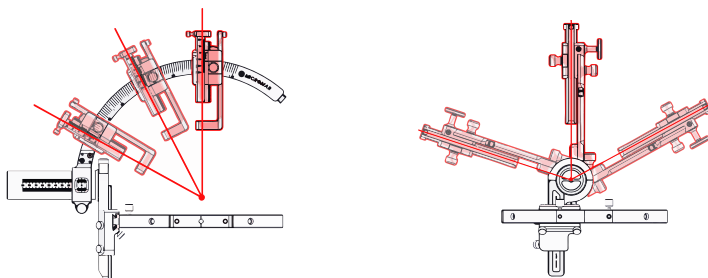
Figure 2

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Technical information on the AimSystem TM03B equipment:

Net weight:

- 2,0kg;

Height:

- 400mm;

Width:

- 422 mm;

Isocenter:

- 243 mm;

Precision:

- Linear precision: 0.5 mm;
- Angular precision: 1°.

8. INTENDED USE

Micromar's Stereotactic AimSystem TM03B makes intracranial accesses feasible using neurosurgical procedures. The Aimsystem TM03B also has a variety of accessories to help the neurosurgeon during surgery in the following situations.

1. In preparing the procedure, stabilizing the patient's position for imaging exams, such as Magnetic Resonance and Computerized Tomography..
2. Targets' locations obtained through the patient's imaging exam.

9. PROCEDURES

The Stereotactic AimSystem TM03B was developed to serve and assist in performing the following procedures:

Guided Neurosurgery:

- Guided surgery;
- Catheter implantation;
- Guided neuroendoscopy;;

Diagnostic Neurosurgery:

- Biopsy;

Therapeutic Neurosurgery:

- Haematoma Evacuation;
- Abscess Evacuation;

Functional Neurosurgery:

- Treatment of movement disorders (Lesioning and Deep electrode-DBS Implantation);
- Treatment of chronic pain
- Behavioral Surgery
- Spasticity
- Treatment of Epilepsy

9.1. COMPONENT LIST

The Figures below show the Basic Components of the Stereotactic AimSystem TM03B.

- **Sterilization Tray 1**

Sterilization Tray 1 is comprised of:

Item	Name	Code	Quantity
1	Sterilization Tray 1	TM0010CAA	1
2	Frame	TM0020FAA	1
3	Reusable Fixation Screws (Stainless Steel)	TM0230FAA	1
4	Reusable Fixation Screws (White Polyacetal)	TM0126FAA	1
5	Extension Posts	TM0140FAA	1
6	Allen Key 4.0mm	TM0032FAA	2
7	Allen Key 6.0mm	TM0034FAA	2
8	Curved Segment	TM0020FAD	1

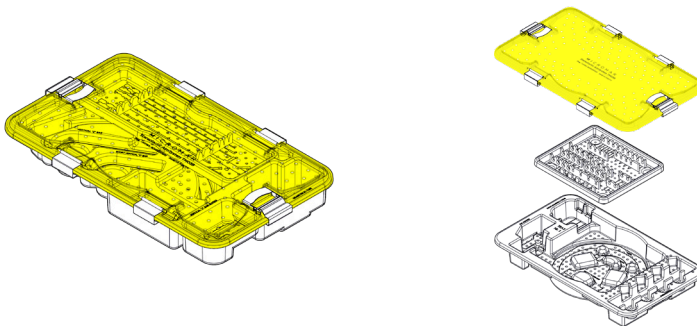
Figure 3

- Sterilization Tray 2**

Sterilization Tray 2 is comprised of:

Item	Name	Code	Quantity
9	Sterilization Tray 2	TM0020CAA	1
10	Millimetrical Ruler	TM0024FAA	1
11	Fixation Screw	PA0020FAA	3
12	"Z" Guide Post	TM0011FAA	1
13	"Z" Post	TM0113FAA	1
14	Short "Z" Guide Post	TM0012FAA	1
15	Short "Z" Post	TM0114FAA	1
16	Millimetrical Tube	TM0122FAA	1
17	Arc	TM0025FAA	1
18	Chariot	TM0010FAA	1

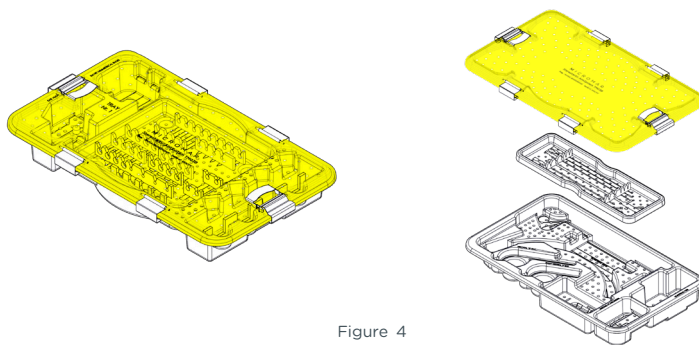
Figure 4

- Support for Basic Computerized Tomography**

Item	Name	Code	Quantity
19	Support for Basic CT	CT0100FAA	1

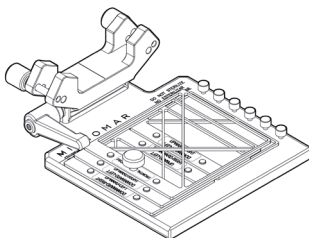
Figure 5

9.2. PREPARING THE PATIENT FOR THE SURGICAL PROCEDURE

To correctly attach the Frame to the patient, the position of the target needs to be verified using a previous imaging exam study.

The Frame should be positioned above or below the target, never in the region of interest, since this will prevent it from being viewed on magnetic resonance and/or computerized tomography exams.

Attachment of the Frame to the patient's head can basically occur in two ways, using DIRECT FIXATION, with the Frame attached to the skull by screws that are screwed directly to the Frame, or INDIRECT FIXATION, using Extension Posts (Optional) that attach the Frame to the patients head using screws..

- **Place the Frame on the patient's head**

Figures (a) and (b) show how accurately the Frame is positioned to perform the surgical procedure.

When the target to be used is located above the Frame, the UPWARD (a) procedure should be followed; and when the target to be used is below the Frame, the DOWNWARD (b) procedure should be followed.

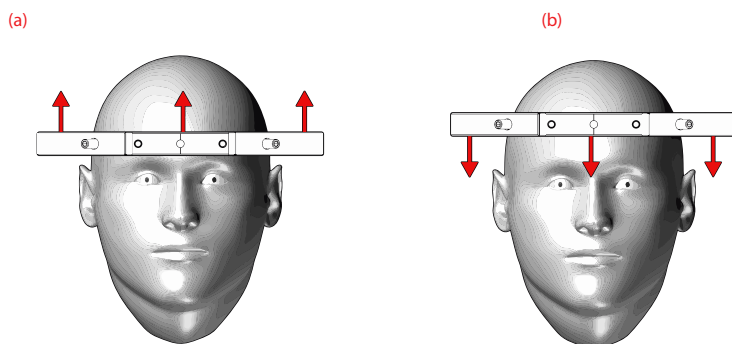
Figure 6

⚠ WARNING: Micromar recommends that the Frame be positioned ± 2 cm above or below the target to be used during the surgical procedure;

⚠ WARNING: The Frame should never be positioned at the same height as the region of interest;

9.2.1. DIRECT FIXATION

DIRECT FIXATION of the Frame on the patient's head can occur using 3 or 4 fixation points, as follows:

- **Direct fixation of the Frame using three screws**

For direct fixation using three screws, the screws are screwed directly into the Frame using the Straight Segment.

If the surgical approach is made using the Ruler attached to the front region of the Frame, a medium or small screw needs to be used on the front part (depending on the size of the patient's head), using two medium or large screws on the back part. Care should be taken so that the front screw does not extend beyond 3 mm from the face (A) of the Frame; otherwise, it will prevent the Ruler from being appropriately placed on the front region.

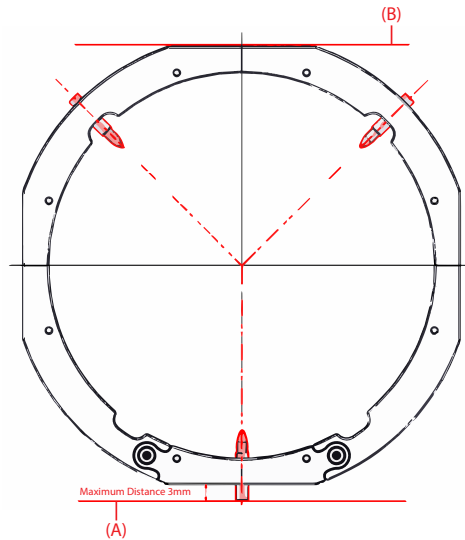
Figure 7

⚠ WARNING: Fixation screws should not extend beyond faces (A) and (B);

⚠ WARNING: Fixation screws have a piercing tip and should therefore be handled with care.

- **Direct fixation of the Frame using four screws**

For direct fixation using four screws, two front screws and two back screws will be used. Screws should not extend beyond faces (A) and (B) of the Frame, since this prevents the Ruler from being attached correctly.

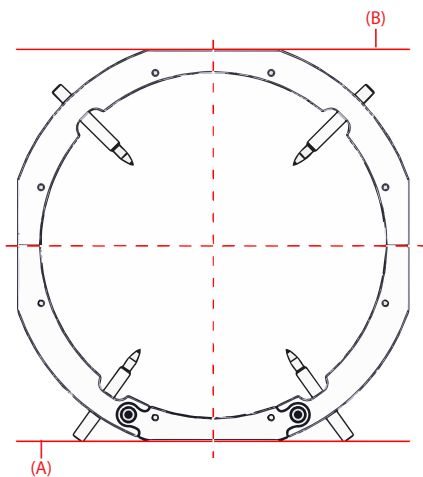
Figure 8

⚠ WARNING: Fixation screws should not surpass faces (A) and (B);

⚠ WARNING: Fixation screws have a piercing tip and should therefore be handled with care.

9.2.2. INDIRECT FIXATION

The Stereotactic Aimsystem TM03B also allows for indirect fixation that allows the Frame to be positioned lower on the skull, which allows targets located at the base of the brain to be reached.

- **Indirect fixation of the Frame using four screws**

For indirect fixation using four screws, the Curved Segment and Extension Posts will be used (in addition to the Frame):

- CURVED SEGMENT

Recommended use

The Curved Segment is recommended when a patient needs to be intubated and the Frame is directly or indirectly attached with four screws. The need to intubate the patient should be assessed in advance, since the system's assembly can not be changed during the surgical procedure.

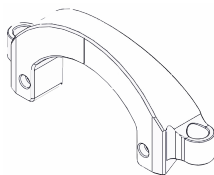
Figure 9

Curved Segment Assembly:

1. To assemble the Curved Segment on the Frame, the Straight Segment must first be removed;
2. After removing the Straight Segment, insert the Curved Segment so that it easily fits;
3. After placing the Curved Segment, tighten the screws (the same ones removed from the Straight Segment), turning them clockwise by using a 4mm Allen Key;
4. To remove the Curved Segment, release the screws holding the segment by using a 4mm Allen Key, turning it counterclockwise.

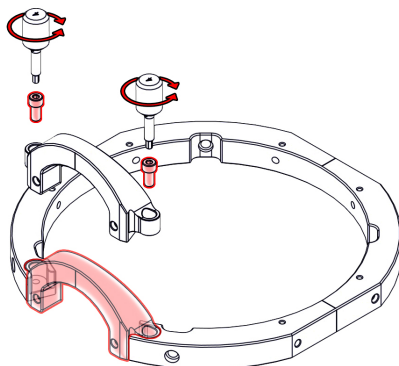
Figure 10

⚠ WARNING: Never change the segment when the Frame is attached to the patient's head, since this will cause irreversible damage the Frame;;

⚠ WARNING: If the Curved Segment needs to be positioned downward, turn the Frame assembled to the Curved Segment before attaching it to the patient's head;

- EXTENSION POSTS

Recommended use

The Extension Posts are used for surgical approaches (always Upward), with indirect fixation, where the target is found on the more basal region of the brain; they are also used in procedures where a greater brain region range is desired.

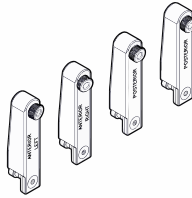
Figure 11

Assembly of the Extension Posts

1. Remove the screws from the Extension Posts, turning them counterclockwise;

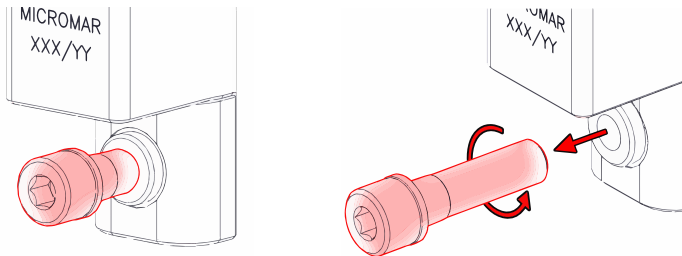
Figure 12

2. Insert the Extension Posts, each in their position, in the Frame's internal part (Figure 13);
3. Lock the Frame posts using the screws (removed in the first step) and use the 4mm Allen Key to attach the Extension Posts to the Frame;

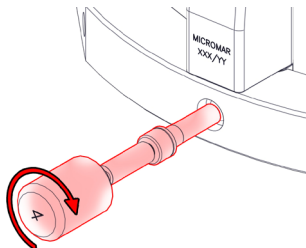
Figure 13

4. If the Reusable Stainless Steel or Disposable Aluminum Screws are used, then the Adapter Bushing Guides need to be used (in this case, use the 4mm Allen Key to tighten the screws);

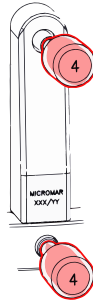
Figure 14

5. To assemble the Adapter Bushing Guides, they must be turned clockwise on the inside of the Frame and the screws should be inserted from the outside to the inside.

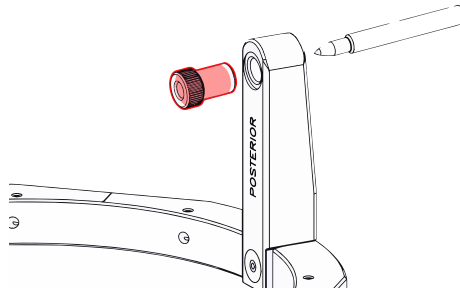
Figure 15

6. If the Reusable Polyacetal Screws are used, remove the threaded Bushing Guides need to be used (in this case, use the 6mm Allen Key to tighten the screws).

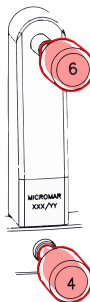
Figure 16

⚠ WARNING: Fixation Screws have a piercing tip and should therefore be handled with care.

- **Indirect fixation of the Frame using three screws**

For indirect fixation using three screws, the A Segment and Extension Posts will be used (in addition to the Frame):

- A SEGMENT (OPTIONAL ITEM)

Recommended use

The "A" Segment is used in conjunction with the Posterior Extension Posts and is recommended when indirect fixation needs to be done using three fixation screws (two posterior and one anterior), allowing for approaches on the frontoparietal region.

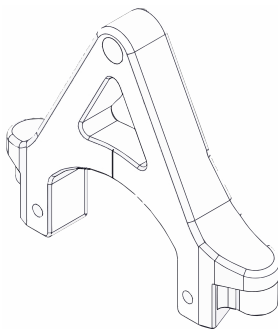
Figure 17

NOTE: The A Segment is an Accessory item; information on use is contained in the Accessories Instruction Manual.

9.2.3. PREPARATION FOR IMAGING EXAM

Micromar has accessories to help the user in obtaining images on equipment for:

- Localizers for Computerized Tomography

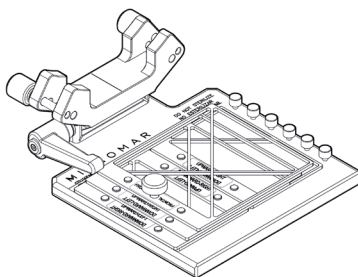
Figure 18

- Localizers for Magnetic Resonance (OPTIONAL ITEM);

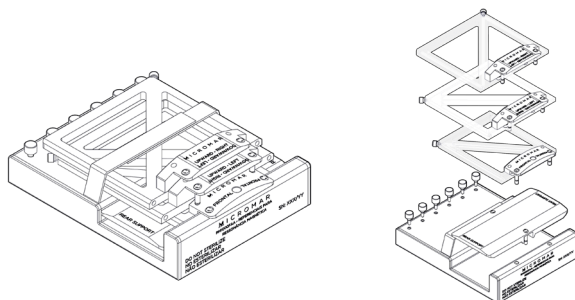
Figure 19

NOTE: The Magnetic Resonance Localizers are Accessory items; information on use is contained in the Accessories Instruction Manual.

9.2.4. MANUAL CALCULATION OF STEREOTACTIC COORDINATES

Stereotactic calculations can be done manually, based on a computerized tomography exam.

⚠ WARNING: Find the exam image that best shows the target or region of interest.

Image example (Upward)

Figure 20

- Determination of the “X” axis**

- With the help of a Ruler, draw one line connecting points 1 and 3 and another connecting points 7 and 9;

Figure 21

2. Calculate and mark the mid-points (A and B) for each of these lines drawn;

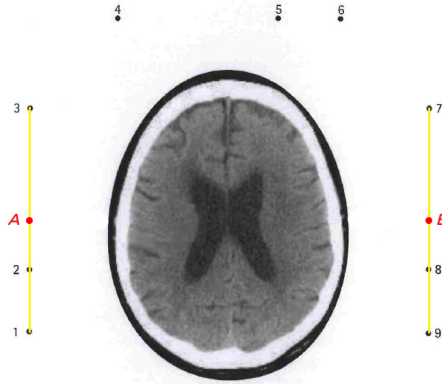
Figure 22

3. Draw a straight line connecting the points (A and B); this line is, therefore, the X axis.

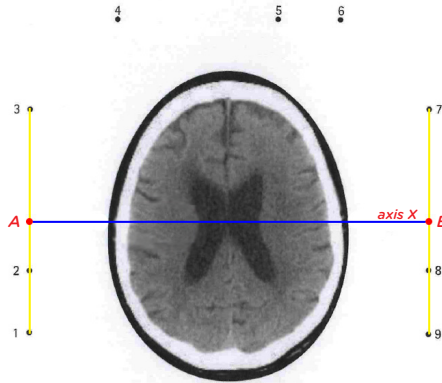
Figure 23

- **Determination of the “Y” axis**

1. With the help of a Ruler, calculate and mark the mid-point (C) of the X Axis (drawn in the previous step);

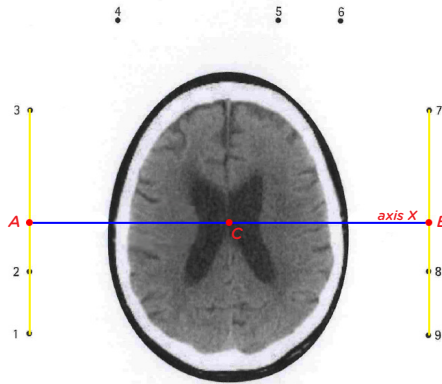
Figure 24

2. Draw a line connecting points 4 and 6;

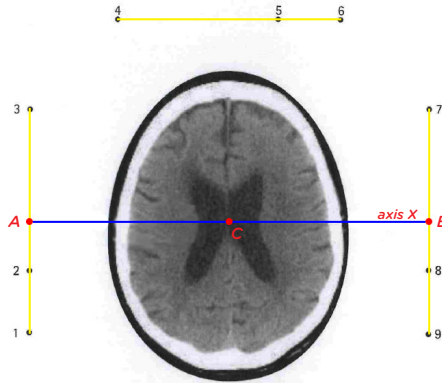
Figure 25

3. Calculate and mark the mid-point (D) of the straight line drawn in the previous step;

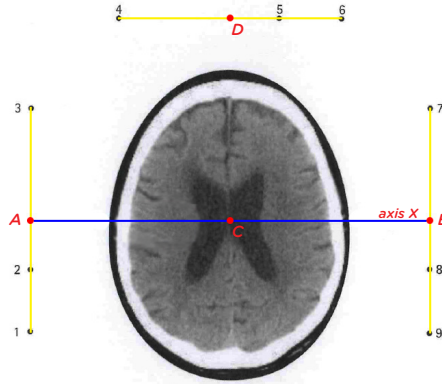
Figure 26

4. Draw a straight line that passes through the mid-point (D) from the previous step and that passes through the mid-point (C) of the X axis. This line therefore represents the Y axis.

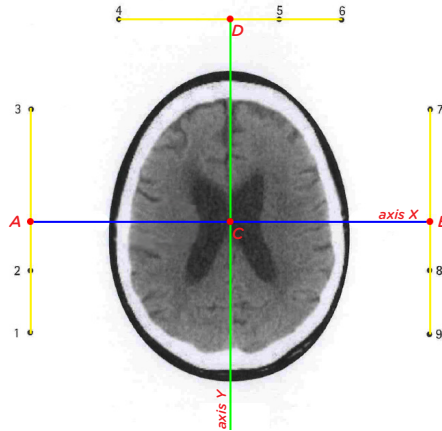
Figure 27

- **Determination of the “Z” coordinate**
- **For “UPWARD” target (positive Z Coordinate)**

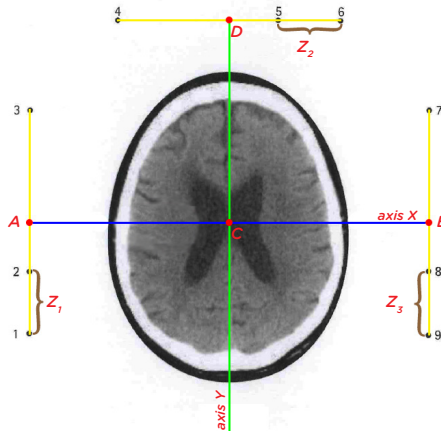
Figure 28

- For targets in the RIGHT POSTERIOR quadrant, measure the distance between points 1 and 2 (Z_1);
- For targets in the LEFT POSTERIOR quadrant, measure the distance between points 8 and 9 (Z_3);
- For targets in the RIGHT and LEFT ANTERIOR quadrant, measure the distance between points 5 and 6 (Z_2);
- **For “DOWNWARD” target (negative Z Coordinate)**

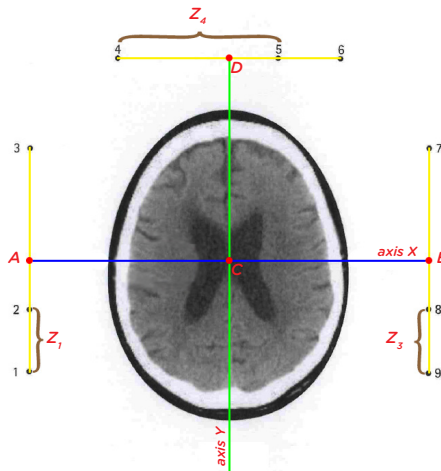
Figure 29

- For targets in the RIGHT POSTERIOR quadrant, measure the distance between points 1 and 2 (Z_1);
- For targets in the LEFT POSTERIOR quadrant, measure the distance between points 8 and 9 (Z_3);
- For targets in the RIGHT and LEFT ANTERIOR quadrant, measure the distance between points 4 and 5 (Z_2);

- **Determination of the target's "X" and "Y" Coordinates**

The intersection of the lines drawn previously is the 0 point in "X" and 0 point in "Y." Therefore, based on these lines, we are able to locate the target coordinates, with the following steps:

1. With the help of a Ruler, draw a straight line parallel to the "X" axis, going through the target; draw a straight line parallel to the "Y" axis, going through the target;

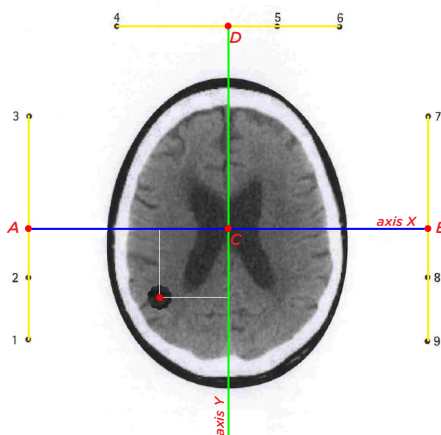
Figure 30

2. The intersection of these two new lines should be the target point (point E);

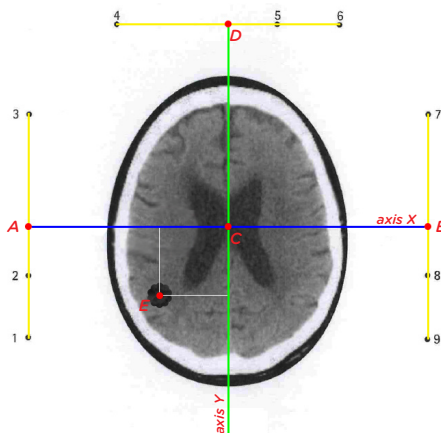
Figure 31

3. Measure the distance (X_1) between point (E) and the "Y" axis. This distance is the "X" coordinate;

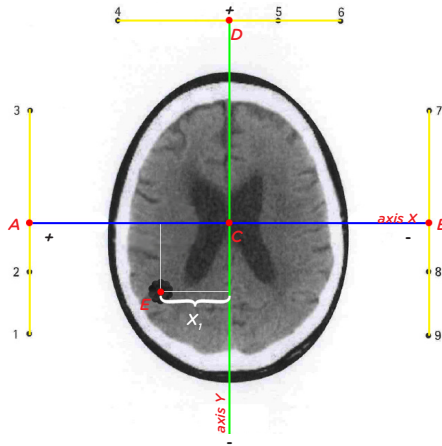
Figure 32

4. Measure the distance (Y_1) between point (E) and the "X" axis. This distance is the "Y" coordinate;

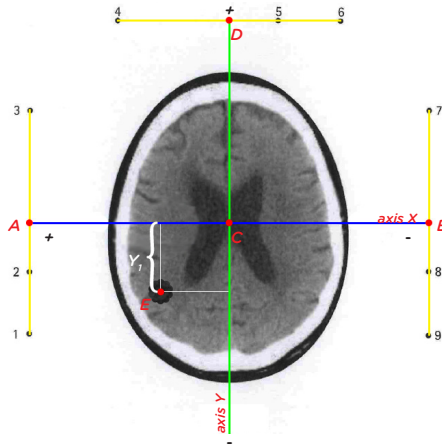
Figure 33

NOTE: These coordinates correspond to the virtual coordinates, which should be converted into real coordinates using the method described in the next topic.

- **Conversion of virtual coordinates into real coordinates**

Print images are not always in real scale (1:1), since the properties can be changed to enter a greater number of images in a single print, causing reduction in these images; with this being the case, dimensions found in the print exam (which can be called virtual coordinates) should be processed and then entered into the equipment.

Virtual coordinates should be corrected using a conversion factor, which is found by dividing the Localization's Real Measurement (containing 140mm) by the Localization's Virtual Measurement, found by using a Ruler to measure the distance between points 1 and 3 in the print exam. As shown below:

$$\frac{140}{\text{MVF}} = \text{CONVERSION FACTOR}$$

After finding the conversion factor, the measurements for the virtual coordinates "X", "Y" and "Z" need to be multiplied by the value found, thus providing the real stereotactic coordinates.

- **Example:**

Suppose that the (virtual) measurements found in the CT exam were "X" = 22mm, "Y" = -18mm and "Z" = 18mm. And the measurement found between point 1 and point 3 was 73mm.

$$\frac{140}{73} = 1.92$$

The result (1.92) found is the conversion factor. Multiply the virtual coordinates found by the calculated conversion factor, according to the table below:

Axis	Virtual coordinate	Conversion factor	Stereotactic coordinate
X	22mm	x 1.92	42,24mm
Y	-18mm		-34,56mm
Z	18mm		34,56mm

The location of the target should be shown by the (+) and (-) minus signs:

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Z Coordinate: Top (+) and Bottom (-).

9.2.5. COMPONENT ASSEMBLY

⚠ WARNING: Prior to initiating assembly of the components of tray #2, the Micromar (HeadFix) fixation system needs to be assembled with the Stereotactic AimSystem TM03B adapter;

⚠ WARNING: For more information on the Micromar fixation system and HeadFix adapter for the Frame, see the equipment user manual for associated use;

⚠ WARNING: Sterilize the Frame before starting the procedure;

⚠ WARNING: All parts should slide together easily; if this is not the case, see the chapter on lubrication (chapter 12);

⚠ WARNING: If there is any resistance, check that the part is correctly assembled or whether a fixation screw is causing interference;

⚠ WARNING: If the problem persists, contact Micromar or an authorized representative in your region;

⚠ WARNING: Fixation screws for each component should be screwed in before following the step by step instructions below.

- **Ruler Assembly**

1. Place the Frame Ruler;

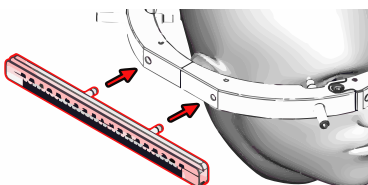
Figure 34

NOTE: The Ruler can be attached to the four Frame quadrants (right and left sides and anterior and posterior sides), establishing its position according to the target location, as chosen by the neurosurgeon.

2. Ensure that the Ruler is lying entirely against the Frame;

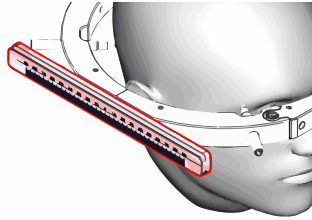
Figure 35

3. Place the two screws and turn clockwise to lock the Ruler in place.

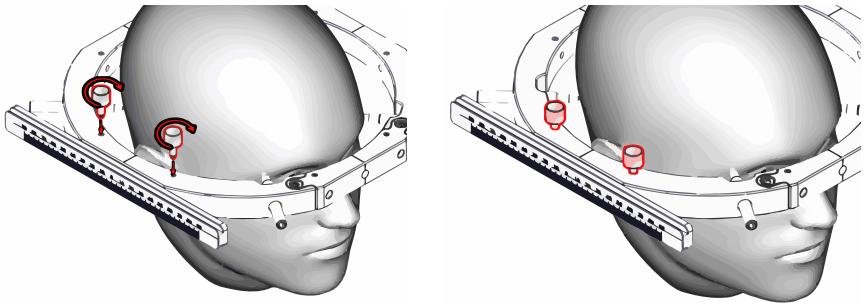
Figure 36

NOTE: When the Upward approach will be used, the fixation screws should be placed in the top part of the Frame, and when the Downward approach will be used, the fixation screws should be placed in the bottom part of the Frame.

- **“Z” Guide Post Assembly**

1. Release all of the fixation screws from the post, so that no interference occurs during assembly;

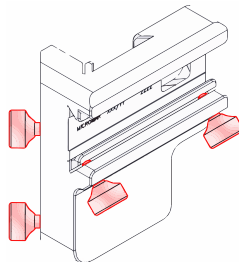
Figure 37

2. After releasing all of the fixation screws, insert the "Z" Guide Post into the Ruler;

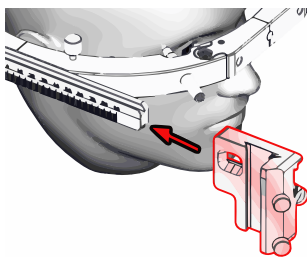
Figure 38

NOTE: When the "Z" coordinate is less than 60mm, the short "Z" Guide Post and the Short "Z" Post can be used.

3. Move the "Z" Guide Post to the desired coordinate;

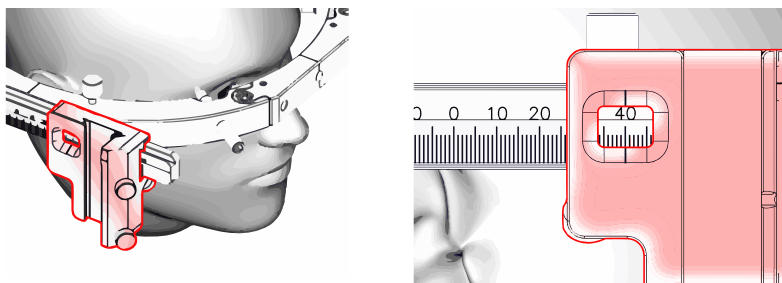
Figure 39

4. After finding the desired coordinate, use the fixation screws to lock the guide in place, turning it clockwise.

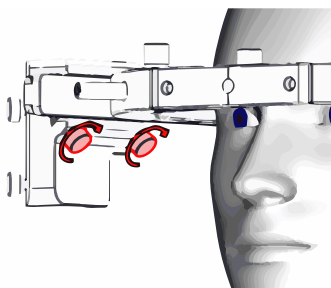
Figure 40

- **“Z” Post Assembly**

1. Release the fixation screw from the “Z” Post, turning it counterclockwise so that there is no interference during assembly;

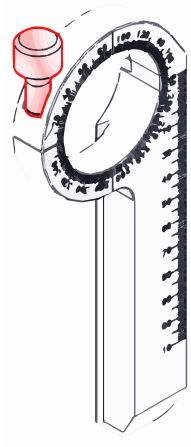
Figure 41

2. After releasing the fixation screw, insert the “Z” Post (the short “Z” Post can also be used with the short “Z” Guide Post, optional items) into the “Z” Guide Post;

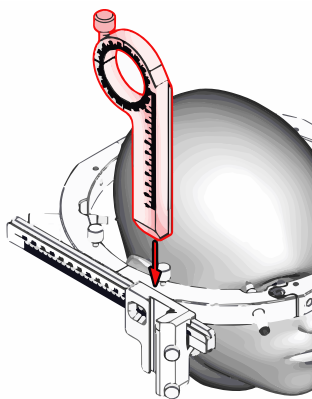
Figure 42

3. Slide the "Z" Post to the desired coordinate;

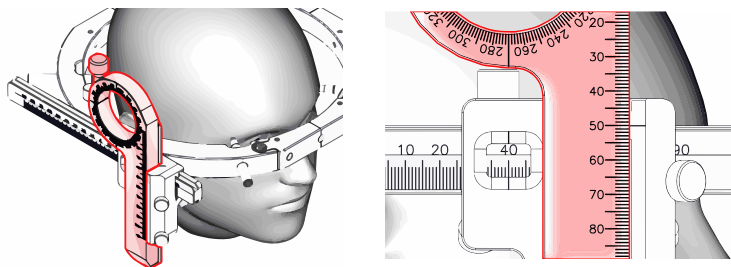
Figure 43

4. After finding the desired coordinate, use the fixation screws on the "Z" Guide Post to lock the "Z" Post in place, turning it clockwise.

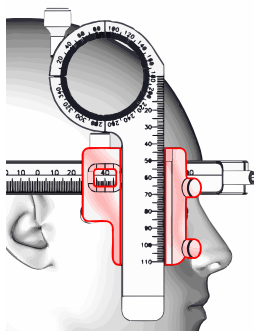
Figure 44

- **Verifying the system's isocenter**

1. Insert the Arc into the Tube;

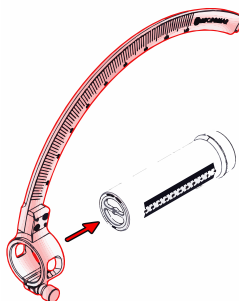
Figure 45

2. Position the Arc on the desired coordinate and lock it into position by using the fixation screws;

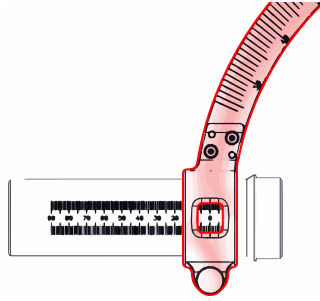
Figure 46

3. Insert the Chariot into the Arc;

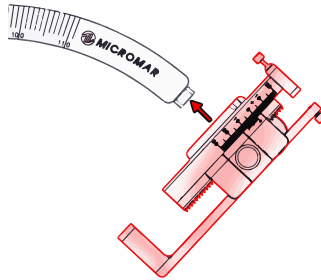
Figure 47

4. Position the Chariot on the "90°" angle coordinate and lock it into position by using the fixation screws;

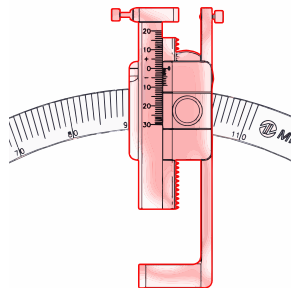
Figure 48

5. On the Chariot's Millimetric Body, go to the "0" coordinate and lock it into position by using the fixation screws;

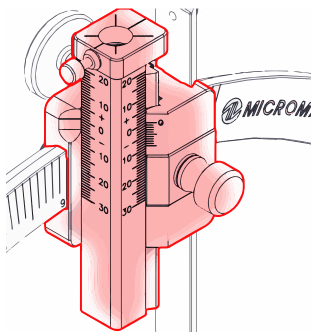
Figure 49

6. Turn toward "LOCK" and lock in place, so that the Chariot runs no risk of falling;

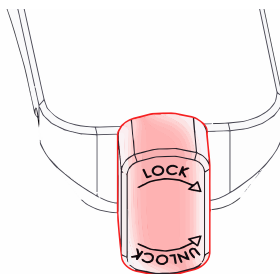
Figure 50

7. Place the Bushing Guides in the Chariot, screwing them in clockwise;

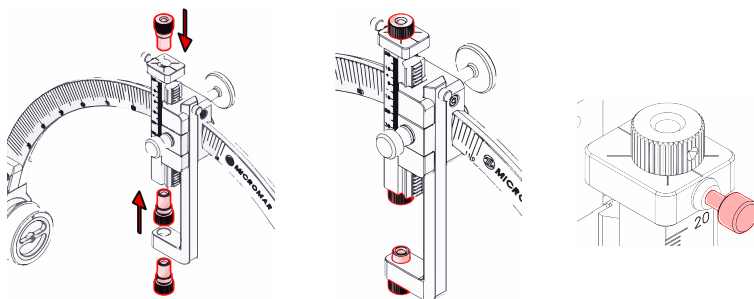
Figure 51

8. Insert the Cannula in the Bushing Guides attached to the Chariot, until it is against the top Bushing Guide;

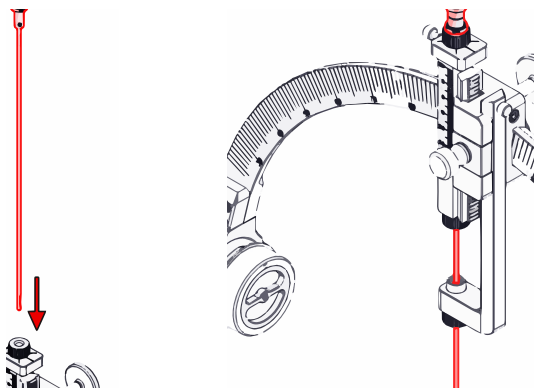
Figure 52

9. Use the Sight on the Tube to check the isocenter according to the cannula model to be used.

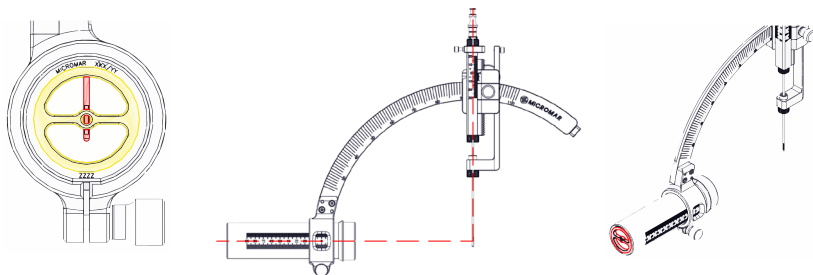
Figure 53

- **Tube Assembly**

1. Remove the Cannula from the Chariot, used to check the isocenter;
2. Check that the fixation screw on the "Z" Post is loose;

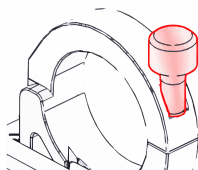
Figure 54

3. Insert the Millimetric Tube in the "Z" Post;

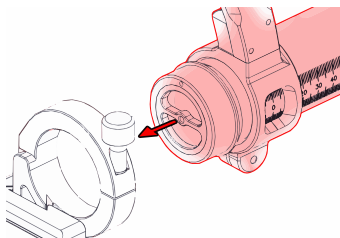
Figure 55

4. Position the Tube in the angular position desired, turning it clockwise or counterclockwise;

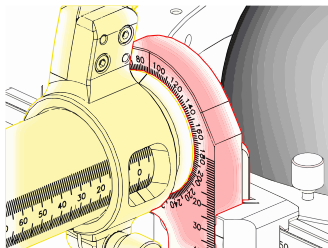
Figure 56

5. Turn the fixation screw clockwise to lock the Tube into place;

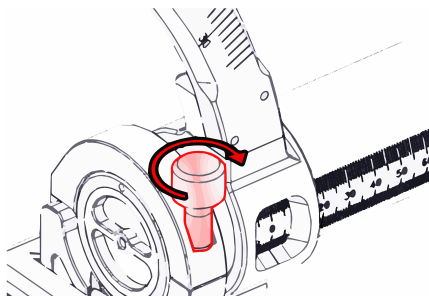
Figure 57

10. SCALES

All of the system's linear scales are graduated in "A" millimeters and all of the system's angular scales are graduated in "B" sexagesimal degrees.

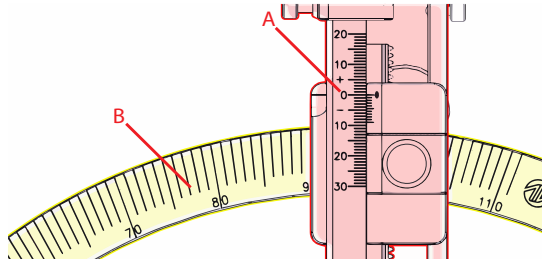
Figure 58

11. TRAYS FOR TRANSPORT / STORAGE / STERILIZATION

To transport the Stereotactic AimSystem TM03B, all parts must be placed in their appropriate locations.

- ⚠ WARNING:** Store in a dry, ventilated and clean environment, protected from directed sunlight and substantial dust;
- ⚠ WARNING:** Three boxes is the maximum amount recommended when piling for transport and storage;
- ⚠ WARNING:** Before using the product, the package should be inspected as to its integrity, in order to guarantee that it is not compromised;
- ⚠ WARNING:** Do not use the material if the package has been damaged and/or violated;
- ⚠ WARNING:** Some of the information found on the package and product label include the maximum number of boxes that can be piled, protect from rain, fragile, and this side up.

12. CLEANING AND STERILIZATION

Cleaning method	
Sanitization	Wet a surgical compress or wipe with enzymatic detergent and apply to the entire surface of the item; rinse using a surgical compress or wipe moistened with water..
Manual Cleaning (immersion)	Immerse the dissembled parts in a container with a solution of cold to warm (up to 40°C) demineralized water and enzymatic detergent (dilute according to manufacturer recommendations). Let the solution act for at least 15 minutes. Using a soft brush, remove residue from the items, paying attention to recesses, screws and lights. Rinse in running demineralized water. Dry using a surgical compress or wipe and clean compressed air.
Automated Cleaning (Ultrasound)	Immerse the dissembled parts in an ultrasonic washing vat with a solution of cold to warm (up to 40°C) demineralized water and enzymatic detergent (dilute according to manufacturer recommendations). Activate the ultrasound for at least 15 minutes. Rinse in running demineralized water. Dry using a surgical compress or wipe and clean compressed air.
Sterilization Method	
VSSP / Autoclave	Sterilization in Saturated Steam under Pressure (SSUP / Autoclave): 134°C ±2 $\frac{2,14 \text{ kgf/cm}^2}{3 \text{ bar}}$ 7 minutes 100% H ₂ O dry for 50 minutes
ETO	Sterilization in Ethylene Oxide (ETO): 50°C ±10 $\frac{-0,3 \text{ bar} \pm 0,1}{-0,3 \text{ kgf/cm}^2 \pm 0,1}$ ≥40%UR 60 minutes 90%ETO and 10%CO ₂ 80 changes/hour ≥ 48 hour of aeration.

The tables below show the cleaning and sterilization methods for the Stereotactic AimSystem TM03B:

Item	Code	Cleaning method			Sterilization Method	
		Sanitization	Manual Cleaning (Immersion)	Automatic Cleaning (Ultrasound)	Saturated Steam under Pressure (SSUP / Autoclave)	Ethylene Oxide (ETO)
Sterilization Tray 1	TM0010CAA	✗	✓	✓	✓	✓
Sterilization Tray 2	TM0020CAA	✗	✓	✓	✓	✓
Support for Basic CT	CT0100FAA	✓	✗	✗	✗	✗
Polyacetal Fix. Screw	TM0126FAA	✓	✓	✓	✗	✓

Key			
✓	Procedure Permitted	✗	Procedure Prohibited

⚠ WARNING: Items exposed to Creutzfeldt-Jakob disease should be discarded according to local regulations.

⚠ WARNING: If the parameters for cleaning and sterilization are not followed as described here, contamination control may not be guaranteed or damages could be caused to the product.

⚠ WARNING: Do not use abrasive and corrosive agents for cleaning, since they can cause damage to the product.

⚠ WARNING: After cleaning, the components should be immediately dried using a dry compress and compressed air. Pay attention to recesses and screws.

NOTE: Follow the washing machine manufacturer's instructions.

NOTE: Follow the information specified by the detergent manufacturer.

13. LUBRICATION

Prior to the surgery, use sterilized vaseline or silicone to lubricate sliding parts and the tips of fixation screws. Lubrication of each part makes tightening and loosening screws easier and facilitates movement of parts during surgery.

14. PREVENTIVE AND CORRECTIVE MAINTENANCE

To guarantee the appropriate functioning and performance of the equipment, it should be returned for preventive maintenance after 18 months.

If there is any defect or malfunction before 18 months have passed, do not use the equipment and immediately contact Micromar or your authorized sales representative.

Micromar is solely authorized to provide technical assistance for the Stereotactic AimSystem TM03B. After the 18 month period has passed, the equipment must be sent every 14 months to Micromar for assessment to find whether any maintenance or disposal is needed.

PREVENTIVE MAINTENANCE TABLE			
Quality Control	Quality Control	Quality Control	Quality Control
1 st Preventive Maintenance	2 nd Preventive Maintenance	3 rd Preventive Maintenance	4 th Preventive Maintenance

15. SYMBOLS

Symbol	Key
	Not sterile
	Keep out of sunlight
	Manufacturer
	Attention, check documents
	Keep out of rain
	Operating manual, operating instructions
	Date of manufacture
	Serial number
	Stacking limit by number
	Fragile
	This side up
	Lot number

16. CARE AND PRECAUTIONS

- The equipment should only be used with the basic component items containing the same serial number (except accessories). Replacement of components with components containing other serial numbers could compromise the equipment's precision, causing damage to the patient.
- Accidental movements involving the Stereotactic AimSystem TM03B components during the procedure can compromise the efficacy of the procedure.
- Avoid dropping, twisting, crushing, impacting or any action that will deform the equipment.
- In the event of any of the above occurrences, the user should contact Micromar for assessment of the equipment.
- Deformation of the equipment's structure compromises the system's precision, placing the patient's life at risk.
- Incorrect reading of equipment scales (user error in reading or because of poor equipment condition) can put the patient's life at risk.
- Screw in and correctly attach the movable parts of the equipment.
- To perform the procedure, the equipment should be cleaned and sterilized according to the instructions in this manual.
- Immediately after the end of the procedure, the equipment should be cleaned according to the instructions in this manual.
- The equipment should be sterilized according to the instructions contained in this User Manual; otherwise, there is a risk of bacterial or viral contamination, causing infectious complications for the patient.
- This product may only be manipulated by a surgeon specialized in neurosurgery and trained by Micromar.
- The surgeon must consider the patient's physical and pathological conditions when attaching the Stereotactic AimSystem TM03B fixation screw to the patient's head, in order to prevent damage to same.
- If the instructions contained in the manual and in the user training are not followed, this could result in incorrect use of the Micromar Stereotactic AimSystem, placing the patient's life at risk.
- Use the maintenance service provided by Micromar or an authorized representative to correctly repair the Micromar Stereotactic AimSystem if any deformations occur.
- To repair or replace any component of this equipment, you must contact Micromar through customer service or through any authorized representative.
- Equipment used to treat patients suffering from Creutzfeldt-Jakob disease MUST NOT BE reused for other patients. Check local disposal laws.
- The Stereotactic AimSystem TM03B is solely meant for use in hospitals and clinics.
- Adverse effects: There are no secondary or side effects directly related to appropriate use of the product.

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QUALITY CERTIFICATION

This manual is aimed at helping the user understand how to correctly, safely and effectively use the Stereotactic Aymsystem TM03B. The user must carefully read the user manual that comes with the product before using it.

This product is warranted by Micromar Ind. e Com. Ltda. within the following conditions:

- 1- An 18 month warranty, as from the date of issue of the receipt of sale to the consumer, for replacement of parts or components that may have manufacturing defects, except for those resulting from improper use or installation, done not following the specifications contained in the user manual;
- 2- Warranty pending presentation of this certificate will be offered at no cost to the customer;
- 3- This warranty is null and void in the event of:
 - 3.1- Defect caused by the customer or by third parties other than the manufacturer;
 - 3.2- The product has been changed or repaired by persons or companies not authorized by the manufacturer;
 - 3.3- The recommendations of the user manual are not followed;

Serial Number: _____

Invoice Number: _____

Date: _____

EQUIPMENT INSPECTED BY:

Quality Control

Product	Code	SN	DatE	PREVENT MAINTENANCE TABLE			
Stereotactic AimSystem TM03B	TM03B	SN: _____	____/____/____	____/____/____ QUALITY CONTROL	____/____/____ QUALITY CONTROL	____/____/____ QUALITY CONTROL	____/____/____ QUALITY CONTROL
				1ST PREVENT MAINTENANCE	2ND PREVENT MAINTENANCE	3RD PREVENT MAINTENANCE	4TH PREVENT MAINTENANCE

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