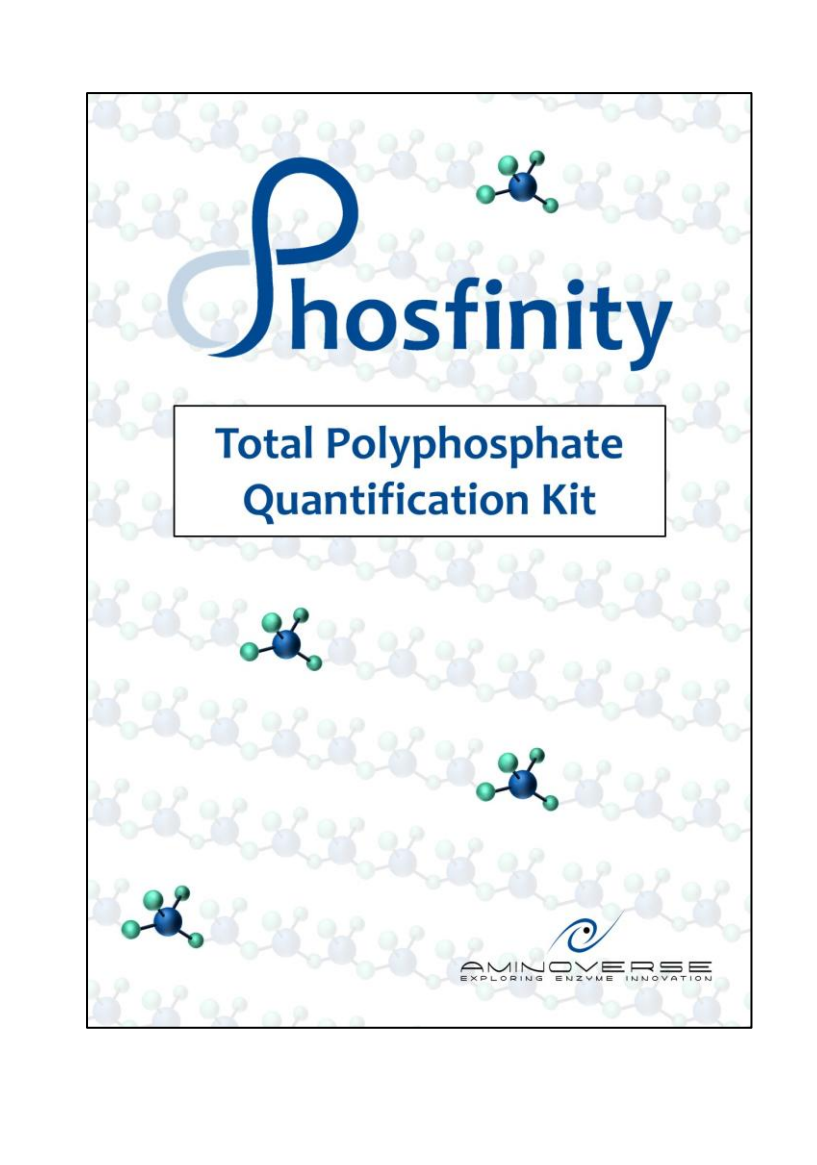


The background of the entire page is a repeating pattern of a molecular structure, likely representing a polyphosphate chain. It consists of a series of light blue spheres (phosphorus) connected by grey rods, with small green spheres (oxygen) attached to the phosphorus atoms. Scattered throughout this background are several 3D ball-and-stick models of a phosphate group, which consists of a central blue sphere (phosphorus) bonded to four green spheres (oxygen) in a tetrahedral arrangement.

Phosfinity

**Total Polyphosphate
Quantification Kit**

The logo for AminoVerse, featuring a stylized blue swirl or 'e' shape.

AMINOVERSE
EXPLORING ENZYME INNOVATION

Contents

1. Introduction	4
2. Key features	5
3. Kit components.....	6
4. Equipment required but not provided	6
5. Storage and stability	7
6. Preparation of reagents.....	7
7. Measuring procedure	10
8. Data analysis.....	13
9. Interfering substances	17
10. Choosing the right extraction protocol.....	18
11. Trouble shooting.....	20
12. Terms and conditions	23
13. References	25

Thank you for your trust and confidence in us and our

Phosfinity: Quant

Total Polyphosphate Quantification Kit.

We always strive to further improve our products, so if you want to give feedback feel free to contact us at:

phosfinity@aminoverse.com.

The kit is sufficient for 100 reactions.

100 samples can be measured if the kit is used in one pass. The kit requires one negative control measurement and one positive control measurement (i. e. 2 reactions) per day.

1. Introduction

Inorganic polyphosphate (polyP, Figure 1) can be found as a linear polymer of phosphate (P_i) in most living organisms¹. Its main features revolve around, e.g., energy storage, blood coagulation, channel formation between cell membranes and biofilm formation^{2,3}. For studies of the phosphate physiology in organisms, it is important to have a reliable and precise quantification method.

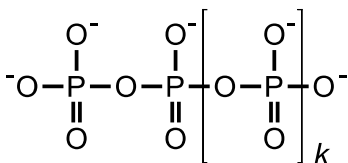


Figure 1: Linear polyphosphate

The new **Aminoverse Total Polyphosphate Quantification Kit** provides a reliable and convenient

option to quantify polyP of all chain lengths, even in highly contaminated samples, down to polyP₂. The assay's high specificity is enabled by enzymatic polyP hydrolysis, followed by colorimetric P_i detection in an easy two step protocol that can be done in less than 1.5 hours with minimal hands-on time.

2. Key features

- High specificity for polyP due to enzymatic detection
- Quantification between 5 to 200 μ M (monomer) polyP
- Detection of all polyP chain lengths down to polyP₂
- Quick colorimetric (2 min) P_i detection. This circumvents artificially high blank measurements due to non-enzymatic hydrolysis of acid-unstable phosphorus containing substances
- Only small amount of sample required (e.g., 200 μ L 100 μ M polyP)
- No interference by commonly encountered substances

3. Kit components

Phosphate standard	5x 1 mL
Negative control	1 mL
Positive control	1 mL
Enzyme	120 µL
Enzyme reaction buffer	10 mL
Phosphate detection solution A	10 mL
Phosphate detection solution B	1 mL

4. Equipment required but not provided

Clear flat bottom 96 well plate

(Multi-channel) Micro pipette

Spectrophotometric multiwell plate reader capable of reading absorbance at 882 nm

5. Storage and stability

The **Enzyme** needs to be stored at **-20 °C** (avoid many freeze/thaw cycles).

Store the remaining kit components at 4-8 °C.

Detection solution A and B should not be exposed to light for extend periods of time.

6. Preparation of reagents

Phosphate detection reagent (prepared freshly)

(50 µL per well, includes negative control, positive control, P_i standards, blanks and samples):

Combine 9 volume parts of solution A and 1 part of solution B (e.g., 45 µL A + 5 µL B). See Table 1 for quick calculation.

Enzyme master mix (always prepared freshly)

(50 μ L per sample. NOT for: P_i standards and blanks):

Combine 49 volume parts of enzyme reaction buffer and 1 part of the enzyme solution (e.g. 49 μ L buffer + 1 μ L enzyme). See Table 2 for quick calculation.

A 20 % excess is recommended to accommodate for loss by using multichannel pipette reservoirs and pipetting error.

Table 1: Preparing the phosphate detection reagent.

Wells	Component [μ]		Wells	Component [μ]	
	A	B		A	B
1	54	6	20	1080	120
2	108	12	30	1620	180
3	162	18	40	2160	240
4	216	24	50	2700	300

5	270	30	60	3240	360
6	324	36	70	3780	420
7	378	42	80	4320	480
8	432	48	90	4860	540
9	486	54	100	5400	600
10	540	60			

Table 2: Preparing the enzyme master mix.

Sample	Buffer [μ L]	Enzyme [μ L]	Sample [μ L]	Buffer [μ L]	Enzyme [μ L]
1	58.8	1.2	20	1176	24
2	117.6	2.4	30	1764	36
3	176.4	3.6	40	2352	48
4	235.2	4.8	50	2940	60
5	294.0	6.0	60	3528	72
6	352.8	7.2	70	4116	84
7	411.6	8.4	80	4704	96
8	470.4	9.6	90	5292	108
9	529.2	10.8	100	5880	120
10	588.0	12.0			

7. Measuring procedure

1. Prepare phosphate detection reagent and enzyme master mix.

📖Note: Store enzyme master mix on ice.

2. Dilute your sample in laboratory-grade water to 5 to 200 μM polyP (as monomer).

📖Note: If the polyP concentration is unknown, test several dilutions.

3. Add 100 μL **negative control** and 100 μL of the **P_i standards** to individual wells on the plate.

📖Note: Use Table 3 as a plate layout example.

4. Add 100 μL of each **sample** into two different wells. Do the same with the **positive control**.

 **Note:** Each sample and positive control requires one well for the actual measurement and one for the blank detecting phosphate contamination.

5. Add 50 μL enzyme reaction buffer (contains no enzyme) to the P_i -standards, one positive control, and the blank measurements.
6. Add 50 μL enzyme master mix to the negative control, the other positive control and the samples.
7. Mix thoroughly by pipetting 3 times up and down (avoid bubble formation).
 **Note:** Set the mixing volume to 90 μL if you have an electric pipette.
8. Incubate at room temperature for 1 h.
9. Add 50 μL of the phosphate detection reagent to every well and mix by pipetting 3 times up and down.
10. Incubate at **room temperature for 2 min.** Remove air bubbles in the wells with fresh pipette tips.

Note: A longer incubation is also possible, but will lead to artificially high blank readings due to non-enzymatic hydrolysis of acid-labile phosphate containing substances (e.g., polyP) in the blank.

11. Read absorption at **882 nm wavelength** in a spectrophotometric plate reader after removing the plate lid.
- 12.
13. Calculate the sample polyP concentration using the Pi calibration curve.

Table 3: Example Pipetting scheme

	Standards	Pi blanks		Negative control	PolyP measurements	
		Positive control	Sample(s)		Positive control	Sample(s)
Well	Well 1-5	Well 6	Well 7 ^a	Well 8	Well 9	Well 10 ^a
Pipetting scheme	100 µl of respective reagent					
	+ 50 µl buffer without enzyme			+ 50 µl buffer with enzyme		
	+ 50 µl Pi detection reagent					
Explanation	Calibration curve	Pi contamination...		Shows the Pi plus polyP	Pi plus polyP concentration...	

		of the positive control	in the sample	contamination of the system	of the positive control	in the sample
--	--	-------------------------------	------------------	--------------------------------	-------------------------------	------------------

^a increase number of wells appropriately for multiple samples.

8. Data analysis

Either follow the instructions below or download the free Excel analysis sheet from:

www.aminoverse.com/Phosfinitivity

Setting up a calibration curve:

Plot the measured absorption values of the P_i -standards (Y-axis) to their known concentration (X-axis; 20, 65, 110, 155, and 200 μM P_i). Draw a linear trendline through the measuring points. Do not include the negative control in

the calibration curve. The formula of the calibration curve (i.e., the trendline) should be in the following form (Eq. 1).

Equation 1:
$$A = S * C + I$$

A: sample absorption

I: trendline intersection with Y-axis

S: trendline slope

C: sample concentration [μM]

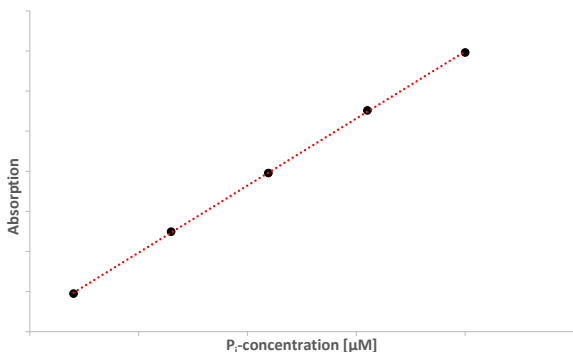


Figure 2: Example standard curve

Calculating the diluted sample concentration

To convert the absorbance values to concentration values, solve the formula of the calibration curve for the variable “C” (Eq. 2). PolyP is detected as its monomer, phosphate, here. The polyP concentration is reported as μM monomer. Calculate the polyP concentration of the positive control, the samples, and the blank measurements. Only values between 5 and 200 μM polyP are valid. Repeat the measurement with higher or lower dilutions if the values are out of the linear range.

Equation 2:

$$C = \frac{A-I}{S}$$

Calculation of the undiluted sample concentration

Subtract the blank value from the corresponding calculated sample polyP concentration and multiply with the sample dilution factor (Equation 3).

Equation 3:

$$\text{true undiluted total polyP} = \\ (\text{sample polyP} - \text{blank sample}) * \text{dilution factor}$$

Quality control

The provided positive control acts as a confirmation of the kit's functionality. The acceptable range for the measured value is 80 – 120 µM polyP.

9. Interfering substances

The following substances do not interfere with the assay:

- ATP up to 2 000 μM
- ADP up to 300 μM
- Glucose 6-phosphate up to 2 000 μM
- Fructose-1,6-diphosphate up to 2 000 μM
- Phosphoenolpyruvate up to 2 000 μM
- NaCl up to 750 000 μM
- NaSO_4 up to 1 000 μM
- ≤ 4.2 g DNA per g polyP (ca. 200 μM DNA monomers)
- ≤ 4.4 g RNA per g polyP (ca. 200 μM RNA monomers)

Contaminating phosphate is detected in the blank measurement. Low polyP concentrations in comparison to the phosphate background level (e.g., < 5 μM polyP per 100 μM phosphate) might be analytically difficult to differentiate from the blank.

10. Using the right extraction protocol

PolyP is often found as an intracellular compound, which needs to be extracted from the cell mass prior to polyP quantification. Even freely dissolved polyP might contain quantification-inhibiting substances (e.g., polyP-degrading enzymes, nucleic acids), which should be removed prior to the polyP analysis.

	Christ protocol⁴	Bru protocol⁵	Werner protocol⁶
Cell lysis	Phenol-chloroform-extraction	Phenol-chloroform-extraction	Inorganic acid (e. g., H ₂ SO ₄)
PolyP purification	None	- DNase/RNase treatment - polyP precipitation with organic solvent	Spin column
Removal of contaminants?	No	Yes	Yes
PolyP recovery	100 %	ca. 60 %	ca. 5 %
PolyP chain length recovered	All chain lengths	Unknown, probably no short chain polyP	> 60 – 80 P-subunits

Suitable for	Cells with > 1 % polyP (as KPO ₃) content in cell dry matter	All sample types	All sample types
---------------------	--	------------------	------------------

Table 4: Comparison of analytical polyP extractions

These three extraction procedures are recommended in combination with the kit (Table 4). However, other extraction procedures are oftentimes similar to the here presented ones and might be suitable as well.

Cells must be lysed to release the polyP. This is achieved with phenol in the Christ and Bru protocols. Phenol must be removed with chloroform later on, but does not degrade polyP. The Werner protocol uses an inorganic acid which can be easily neutralized in higher throughput. Inorganic acid, however, does not lyse the cells fully and hydrolyses polyP.

After polyP release from the cell, the Christ protocol does not further purify the polyP. Therefore, no polyP is lost. In the Bru protocol, nucleic acids are degraded with DNase and RNase and the polyP is precipitated with an organic solvent. Ca. 40 %, including short-chain polyP, of the polyP is lost during the polyP purification. In the Werner protocol, the polyP is bound to a spin column membrane. In this purification, ca. 95 % of the polyP, including polyP with a chain length shorter than 60 to 80 P-subunits, is lost. One advantage of the Bru and Werner

protocols is the high polyP purity that is achieved. PolyP-degrading enzymes and other contaminants, are efficiently removed. The low polyP recovery of the Werner protocol is compensated for by a convenient and quick to perform lab protocol. Concluding, there is a negative correlation between the polyP recovery and polyP purity for polyP extraction protocols. In doubt, the Bru protocol should be tested first. If you are expecting little contaminants or a cellular polyP content of >1 % polyP (as KPO_3) in dry matter, the Christ protocol might be more appropriate because it recovers more polyP.

11. Trouble shooting

Q.: *"The polyP concentration I measure is way below the expected level."*

A.: PolyP-degrading enzymes in your sample or harsh conditions (high temperature, non-neutral pH, higher valent cations, such as magnesium and calcium) might cause polyP degradation. Try a different sample extraction protocols (see the chapter above).

Is your positive control in the range of 80 to 120 μM polyP? If not, there might be an error in assay execution.

Despite its less acidic buffer system compared to competing products, the phosphate detection reagent is still able to hydrolyse acid-sensitive phosphate compounds, such as polyP, after extended periods of time. This can cause an increased blank value and will consequently lead to a falsely low estimation of the sample polyP concentration. Adhere to the 2 min incubation time after phosphate detection reagent addition before reading the absorption.

Q.: *“Is the kit suitable for mammalian cells?”*

A.: Yes. The kit measures the total concentration of dissolved polyP. Therefore, it does not matter from where the polyP stems. Two interfering substances might be faced: A) If the sample contains too much phosphate, the blank reading will be high. If the polyP concentration is low, it might be hard to differentiate the polyP from the background phosphate. B) If you have more than 4.2 g of nucleic acid per g polyP, the nucleic acid will interfere with the measurement. Use an appropriate extraction protocol in both cases (see the chapter above).

Q.: *“I want to measure samples with a predicted high polyphosphate concentration. Why is the highest provided standard at 200 μ M?”*

A.: The detection reagent has a linear range between 5 and 200 μM polyP. Every sample with a concentration above this value needs to be diluted prior to testing.

Q.: *"My sample contains a lot of nucleic acid. Can I use the kit?"*

A.: Up to 4.2 g DNA per g polyP or 4.4 g RNA per g polyP (ca. 200 μM nucleic acid monomer) are tolerated by the assay. Higher nucleic acid concentrations must be removed by an appropriate extraction procedure (see the chapter above).

Q.: *"Can I use freeze-thawed / dried cells?"*

A.: Yes. However, the polyP needs to be extracted from the cells (see the chapter above).

Q.: *"Will the test measure both free and bound polyP?"*

A.: The kit measures all water-soluble polyP that is accessible at least from one end of the polyP chain. Covalent attachment of one polyP end is no problem.

Q.: *"I measure a few samples every day. Why isn't there enough enzyme for the stated 100 samples?"*

A.: The kit procedure includes two controls that also require enzyme. Measuring only a few samples at a time will increase the amount of enzyme used for the controls.

12. Terms and conditions

The Phosfinity kit and its components are intended for research use only, except anything else is stated in an official Aminoverse B.V. document. Furthermore, the Phosfinity kit is suitable for *in vitro* use only and has to be handled by qualified personnel.

The purchase of this product only authorises the buyer to application of the kit for internal research and does not include the right to resell, repackage or sublicense it. The buyer may also refrain from disassembling, modifying, enhancing, reverse engineering,

supplementing or changing the Phosfinit kit or its components, unless stated otherwise. Furthermore, any type of clinical diagnostic or application on the human body is prohibited.

Aminoverse B.V. only warrants the conformity of its products to the stated specifications and performance ranges. Aminoverse B.V. is not liable for any damages resulting from improper applications, test environments or storage conditions. Furthermore, no claims for any damages, whether direct or indirect, incidental, compensatory, foreseeable, consequential, whether based upon warranty, contract or strict liability arising in connection with the sale or the failure of the Phosfinit kit or its components to perform in accordance with the stated specifications.

Aminoverse B.V.'s sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event the kit fails to perform as warranted.

Aminoverse B.V. is not liable for damages or defects arising in shipping and handling, or out of accident or improper or abnormal use of this product; defects in products or components not manufactured by Aminoverse B.V., or damages resulting from such non-Aminoverse B.V. components or products.

Aminoverse B.V. makes no other warranty of any kind whatsoever, and specifically disclaims and excludes all other warranties of any kind or nature whatsoever, directly or indirectly, express or implied, including, without limitation, as to the suitability, reproductivity, durability, fitness for a particular purpose or use, merchantability, condition, or any other matter with respect to Aminoverse B.V. products.

©Aminoverse B.V. All rights reserved. All patents and trademarks are the property of Aminoverse B.V.

13. References

1. Gerasimaite R, Mayer A. Enzymes of yeast polyphosphate metabolism: Structure,enzymology and biological roles. *Biochem Soc Trans.* 2016;44:234-239. doi:10.1042/BST20150213

2. Munees A, Almas Z. Biological Importance of Phosphorus and Phosphate Solubilizing Microbes. 2009;(April):1-14.
10.3389/fmicb.2020.588605
3. Takeda E, Taketani Y, Sawada N, Sato T, Yamamoto H. The regulation and function of phosphate in the human body. *BioFactors*. 2004;21(1-4):345-355. doi:10.1002/biof.552210167
4. Christ JJ, Blank LM. Analytical polyphosphate extraction from *Saccharomyces cerevisiae*. *Anal Biochem*. 2018;563:71-78. doi:10.1016/j.ab.2018.09.021
5. Bru S, Jiménez J, Canadell D, Ariño J, Clotet J. Improvement of biochemical methods of polyP quantification. *Microb Cell*. 2017;4(1):6-15. doi:10.15698/mic2017.01.551
6. Werner TP, Amrhein N, Freimoser FM. Novel method for the quantification of inorganic polyphosphate (iPoP) in *Saccharomyces cerevisiae* shows dependence of iPoP content on the growth phase. *Arch Microbiol*. 2005;184(2):129-136. doi:10.1007/s00203-005-0031-2

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

This image shows a single page of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

Aminoverse B.V.

Daelderweg 9

Nuth, 6361HK

Tel.: +31 452 084815

Mail: phosfinity@aminoverse.com

www.aminoverse.com